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Athlete mode

☐

The personal information is only assessed as health data, Fitdays will keep you confidential

Confirm

Male

Female

Nickname

Fitdays

Birthday

29/11/1995

Height

163 cm

Target weight

54.7kg

Athlete mode

☐

The personal information is only assessed as health data, Fitdays will keep you confidential

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GUARANTEE DOCUMENT

Thank you for purchasing this COMELEC product, this device has been checked in all its components meticulously, so we guarantee the correct operation, but for this to be the case we need that YOU read the instruction book carefully, as this will contribute to the best performance of this article.

P.G.E. S.L guarantees to the buyer of this product, a total guarantee of TWO YEARS, governed by the established conditions whose regulation is found in Royal Legislative Decree 1/2007 of 16 December, which approves the revised text of the General Law for the Defense of Consumers and Users, (articles 114 to 127)

P.G.E. S.L is committed to repairs are completely free, labor and spare parts, including.

The Technical Assistance Service will resolve any repairs that may be required by this apparatus.

Carefully review the INSTRUCTION BOOK included, on the use, maintenance and conservation of the device.

Any anomaly that could occur due to the non-observance of such instructions will not be protected for the guarantee.

The date of sale of the device, for the purposes of warranty period, will be that which appears on the purchase invoice.

Guarantee operation: As regards all those compliance defects detected by the technical assistance, which are manifested within six months of the delivery of the goods, it is assumed that they already existed on that date, the intervention will be completely free, (unless this hypothesis is incompatible with the nature or nature of the conformity defect).

During the following eighteen months, it must be the user who demonstrates. That the lack of conformity it was coming at the time of delivery.

MANUFACTURED CHINA FOR COMELEC

DECLARATION OF CONFORMITY

Manufacturer Name: PRAGMATICA GRUPO EUROPEO, S.L

Address: Tapissers 10 Poligono Industrial Els Mollons

Post Code, 46970 City: Alaquas VALENCIA

We declare under our responsibility that the product:

Product Description: Smart body fat scaler

Model ; BBCH6103

Brand: COMELEC

It complies with the European Union legislation in compliance with the 2004/108 / EU guidelines

The standards apply:

EN 61000-3-3:2013

Name: Pablo Gonzalez

Position : Administrador

Date: 10 Enero 2.019

Signature:

INFORMATION ON RECYCLING

Disposal of the old appliance, the European Directive 2002/96 / EC on waste equipment Electrical and Electronic (WEEE), establishes that old appliances can not be thrown as municipal waste without normal classification. Appliances have to be collected separately with the to optimize the recovery and recycling of the materials they contain and reduce the impact on human health and the environment.

The symbol "container with wheels" crossed out on the product reminds you of obligation, that when you dispose of the device must be collected separately. Consumers should contact the local or retail authority to obtain information about the correct disposal of your old device.

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DATE OF SALE

Seller's signature and stamp

IMPORTANT: To make the guarantee effective, the SELLER must fill out, sign and stamp this Certificate. To receive assistance under the conditions set forth in this certificate, please submit to the Technical Service Authorized, this document next to the Purchase Invoice.

TROUBLESHOOTING GUIDE

1. Bluetooth not connected: The scale and smart phone must be both have Bluetooth on. Check the IOS version of your smart phone; it must be IOS 8.0 or higher / Andorid 5.0 or higher. Click Me --> device find the scale. **2. No body fat data measured when weighing.** User must be barefoot to measure the body fat. Make sure your feet and the scale are both dry. Scale not finishing weighing when the number display is flashing (Keep standing on the scale around 10 second until the number on the display stops flashing). Bluetooth is not turned on. Bluetooth is paired with a difference scale. **3. The scale is producing inconsistent weights:** Check whether the scale is on a hard, flat surface. Check each sensor foot on the back to make sure nothing is stuck to the bottom of it. Remove the batteries to restart the scale. **4. The App is display the incorrect flat-muscle ratio.** The body parameter setting is incorrect, check to make sure that you enter the correct gender, height, and age. **5. How do I recalibrate the scale after moving it or change the batteries?.** Step on the scale to turn it on. Let is automatically drop to 0.0kg/lb to calibrate. **6. The scale does not turn on.** Remove the batteries and put it back. Make sure to insert the batteries in the correct direction. **7. How to change unit?.** You need to change the scale unit by the backside button. Then APP will show the same unit with the scale.

BATTERY ADVICE

If the scale is not to be used for long periods, it's recommended to remove the battery to avoid damage to the scale due to possible battery leakage. Keep the battery out of reach of children.

Do not mix old and new batteries, with different compositions or of different brands in order to prevent possible leakage, explosion.

Do not heat or deform the batteries or explore to fire.

Waste batteries should not disposed of with household waste.

Please check with your local authority for battery recycling advice.

Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local Authority or retailer for recycling advice.

DISPOSAL

The packaging material is recyclable. Please dispose of it in an environmentally friendly manner and take it to a recycling centre.



Dispose of the device in an environmentally friendly manner. Dispose of it at a recycling centre for used electrical and electronic devices.

You can obtain more information from your local authorities.

The batteries should be removed prior to disposal of the device and disposed of separately from the device. To protect the environment, non-rechargeable and rechargeable batteries may not be disposed of with normal household waste, but must be taken to suitable collection points. Also note the relevant statutory regulations governing the disposal of batteries.

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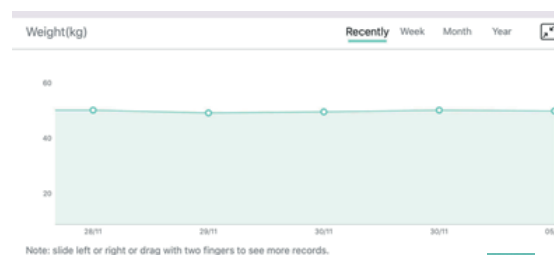
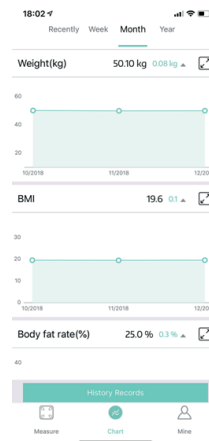
9. Add, switch user by click,  or click Me --> Manage users, then add personal data. The Smart scale can recognize which user is step on it and remind user to switch the correct user. Each account can support 24 users.

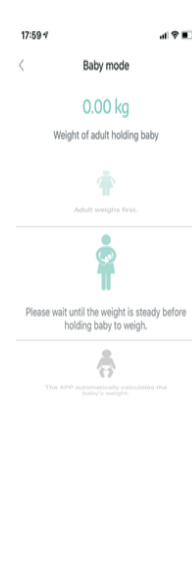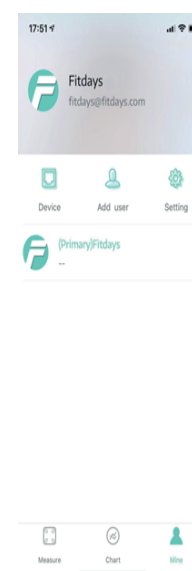10. Weighing baby: click + --> Baby Mode. Step on the scale by your own, after scale gets your weight, and then hold the baby to get the weight of baby.

5. Tap your foot on the scale to activate the display. Wait until the display shows “0.0”(kg/lb)
6. Step on the scale with BARE feet to turn the scale on.

Clean and Bare feet


With socks or shoes

NOTE:

If stepping on the scale with socks, only weight and BMI will be measured and displayed on the main page. To get complete data, please step on the scale with bare feet. 7. The number on the display is flashing when the scale is weighing and analyzing. Continue to stand on the scale until the circle is stable and shows, your phone has received all the data from the scale.

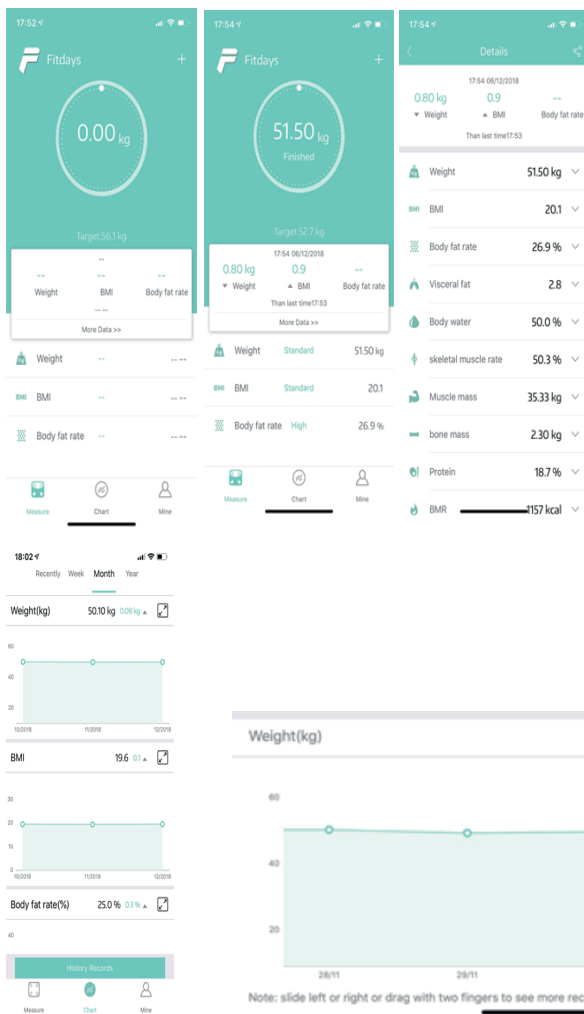8. You can find the analyzing report on Body Weight, BMI, BMR, Muscle, Bone Mass etc. (see full details on More data >>). When you click Record, you can see your long-term track chart form. Fitdays supports you to Zoom + the record for every functions via the right icon; Fitdays also supports record for recently, week, month and year;

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





2. Turn on the Bluetooth on your smart phone or tablet.
3. Resister your own ID by email or social media account.

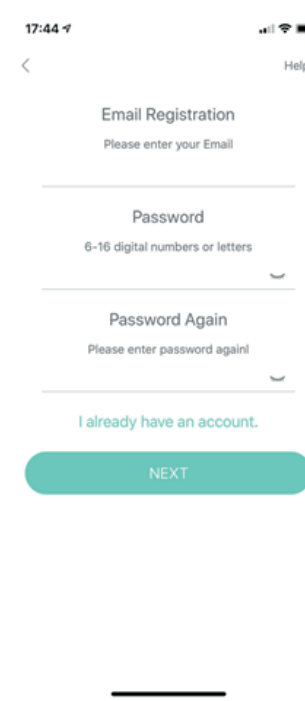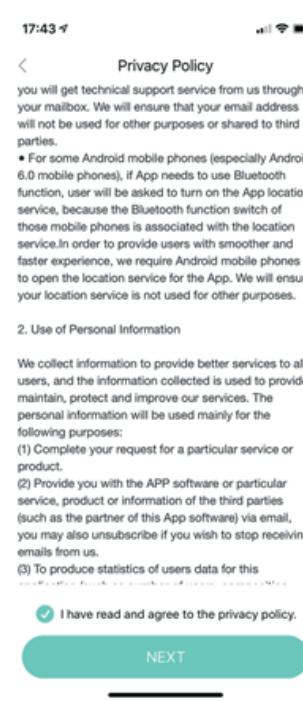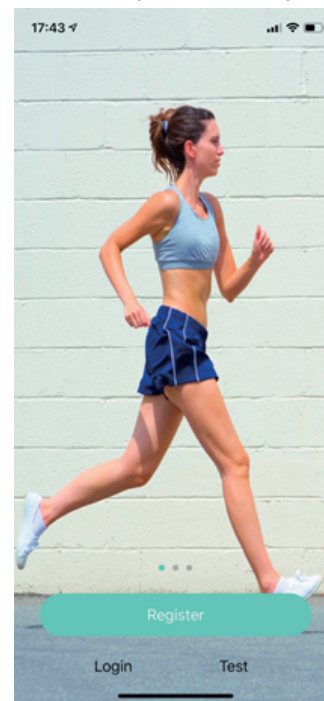4. Add personal data and confirm OK. **Note:** Enter the wrong gender, age, height, or corporeity level may result in measurements being inaccurate.

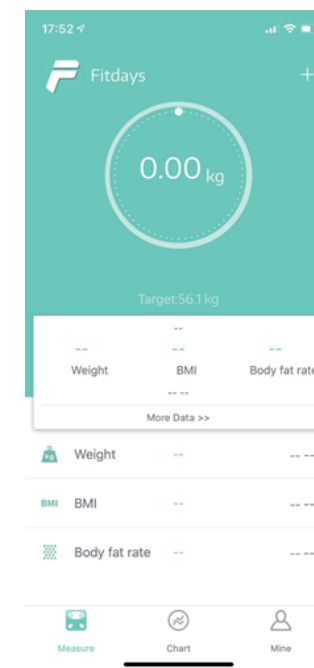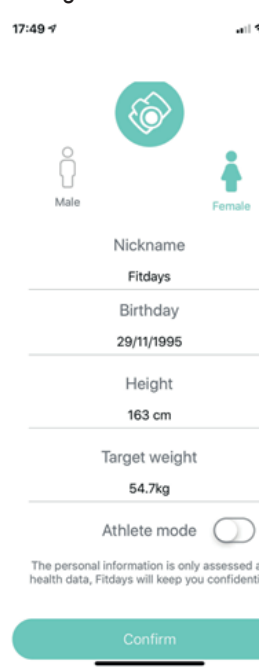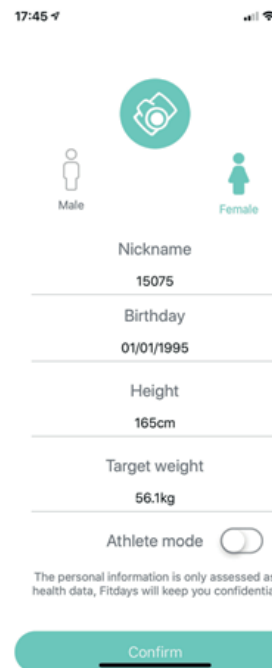Thanks for choosing smart body fat scale. This fully body composition scale can help you conveniently measure body weight, BMI, body fat rate, visceral fat, body water, skeletal muscle rate, muscle mass, bone mass, protein, BMR and body age.

RECORDED PARAMETERS: ●Weight ●BMI ●Body fat rate ●Visceral fat ●Body water ●Skeletal Muscle rate ●Muscle mass ●Bone mass ●Protein●BMR ●Body age ●Height record (<16 years old)

INDICATION SYMBOLS

Lo	Low battery – please change the batteries
Err	Overload
C	An error has occurred during measurement. Step off the scales and step on them again to repeat the measuring process.

WHAT IS IN THE PACKAGE: 1 x Smart Scale 1 x User Manual

WARNING

DO NOT use with medical implants such as pacemakers. **DO NOT** stand on the edge of the scale or jump on it. **DO NOT** overload the scale (MAX. 180kg / 396lb). **DO NOT** drop scale or droop and objects on it as this may damage the sensors. **DO NOT** immerse the scale in water or use chemical cleaning agents. Clean the scale with slightly damp cloth. **NOT** recommended for pregnant women. **NOT** recommended for infants, toddlers, and children under 10 years of age. Any measurement obtained using this device is for reference only and should not be considered as a medical opinion. Consult your primary care doctor or physician before making changes to your diet, exercise plan or physical activities. Always place the scale on a hard, dry and flat surface before measurement. Make sure your feet are dry before stepping on the scale. Connect the batteries in accordance with the correct polarities indicated. Remove batteries if the scale is not be used for a prolonged period of time. The user under 6 years old, only weight function, no BMI and body fat, etc. data. The user under 10 years old, only weight & BMI function, no body fat, etc. data. **This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. children shall not play with the appliance. cleaning and user maintenance shall not be made by children without supervision.**

SPECIFICATION

LCD display:

Weight unit: kg/lb (switch the units in the bottom of the scale body)
Weight Limit: 6kg – 180kg (13.2lb – 396lb)
Account number: No limitation (each account can add 24 persons)
Baby mode: Allow you to scale the weight of baby
SUPPORT DEVICES: IOS 8.0 or Higher Android 5.0 or Higher
INSTALL FITDAYS: 1. Search “Fitdays” in Apple Store or Google Play or scan the QR code below: 2. Download and install the App on your device. 3. The App icon will appear on your phone or tablet after the installation is completed.

HOW TO USE APP: 1. Install the batteries from the package or pull out the insulator if batteries have been inside already. Place the scale on the hard and flat ground.

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INSTRUCTIONS MANUAL SMART BODY FAT SCALE MODEL BBCH 6103

CE

SUITABLE FOR FAMILY USE
THIS PRODUCT IS FOR HOME USE ONLY

PLEASE READ THE INSTRUCTIONS CAREFULLY OF THIS MANUAL BEFORE STARTING.

DO NOT HANDLE THIS SMART BODY FAT SCALE INTERNALLY.

The texts, photos, colors, figures and data correspond to the current technical level from the moment they have been printed. We reserve the right to modifications motivated by the permanent development of the technique in our

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